



# Accessories



## Diagnostic Instruments

### DIAGNOSTIC KIT

CAT#	DESCRIPTION
9-7005-00	<b>Diagnostic Kit includes:</b> Battery handle Otoscope chamber Ear tips: size 2.5, 3.5, & 4.5mm Nasal speculum screw type Ophthalmoscope head Metal tongue depressor Illuminated tongue holder Laryngeal mirror no. 3 Laryngeal mirror no. 4 Tongue holder for wooden blades

CAT#	DESCRIPTION	BOX
5-0240-59	Replacement Lamp for Otoscope Handle	1



### OTOSCOPE

- Standard fittings
- Chrome plated brass handle
- Ear tips are reusable and autoclavable
- Packed in a zippered pouch
- Adapter is available which allows use of Welch Allyn disposable specula tips
- Requires two "C" batteries (not included)

CAT#	DESCRIPTION
9-7004-20	<b>Otoscope, includes:</b> Battery handle Otoscope chamber Ear tips: size 2.5, 3.5, & 4.5mm
9-7004-05	Adapter for Welch Allyn disposable specula tips





# Accessories




## Nerve Stimulators

These units, though small, still offer sufficient output to ensure current is fully adjustable from 0 to 70 mA. Instruction manual, ball electrodes, battery and leads with alligator clips are included with each unit.

### MICRO STIM

- Variable output: 0-70 mA
- **Switch functions:**
  - Twitch: 1 pulse per second
  - Tetanus: 100 Hz
  - Built-in train of four: 4 pulses per 2 seconds
  - Stand-by position: low battery drain
- Pulse indicator light: flashes red when battery needs to be replaced
- Standard battery: 9 volt alkaline
- Lightweight: 7 oz. including battery

CAT#	DESCRIPTION	BOX
8-1053-60	Micro Stim Nerve Stimulator	1

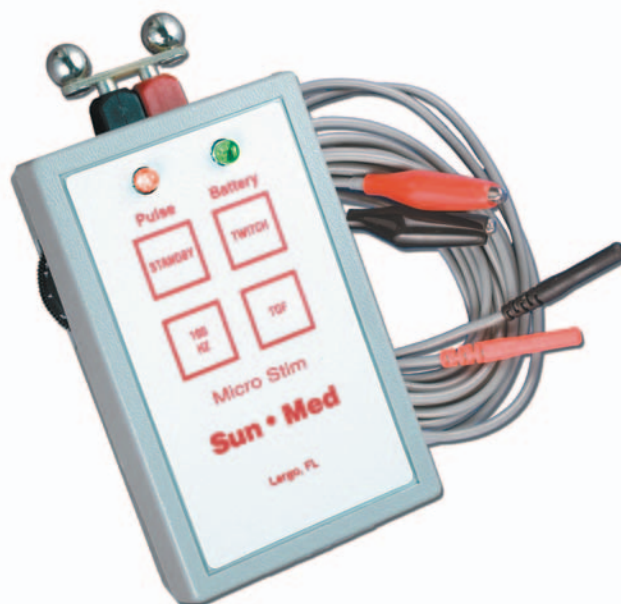

### MICROSTIM PLUS

- Double burst: 2 60-ms bursts of 50 Hz separated by .75 s
- Twitch: 1 s (1 pulse/s)
- Tetanus: 50 Hz (50 pulses/s) or 100 Hz (100 pulses/s)
- Built-in train of four: single (4 pulses/2 s)
- Output current: adjustable 0-70 mA
- Stimulus pulse: square wave monophasic pulses (200- $\mu$ s duration)
- Indicators: battery LED (green) indicates that the unit is on, flashes on low battery
- Pulse LED (red) flashes and audible indicator sounds each time a stimulus pulse is generated.

CAT#	DESCRIPTION	BOX
8-1053-62	Microstim Plus Nerve Stimulator	1



### ACCESSORIES FOR BOTH UNITS

CAT#	DESCRIPTION	BOX
8-1053-76	Stim Ball Electrode & Hollow Pin	Pr
8-1051-78	6' Snap & Shield Pin	Pr
8-1051-79	Stim Lead 6' Clip & Shield Pin	Pr
8-1050-09	Carrying Case	1



# Accessories



## Stethoscopes



### DOUBLE HEAD

- Molded one piece Y tubing
- High sensitivity - lightweight
- Non-chill ring
- Anodized aluminum chestpiece
- Chrome plated brass binaurals
- White plastic eartips
- Superior acoustics
- 31" overall length
- Latex free

CAT#	DESCRIPTION	BOX
6-0027-40	Double Head, Black	1
6-0027-41	Double Head, Red	1
6-0027-42	Double Head, Gray	1
6-0027-43	Double Head, Green	1
6-0027-60	Double Head, Black	10
6-0027-61	Double Head, Red	10
6-0027-62	Double Head, Gray	10
6-0027-63	Double Head, Green	10

### SPRAGUE RAPPAPORT

Eliminates the need for different stethoscopes for infants, children and adults. Selection of the proper head permits desired acoustical responses in low, medium and high pitched sounds. Heavy wall 22" tubing and choice of ear tips help reduce leakage of sound and extraneous noise. This latex free stethoscope comes complete with an accessory kit. Overall length 29".

#### ACCESSORY KIT CONTAINS:

- 1 pair black plastic ear tips
- 3 open bells: infant, child & adult
- 1 pair soft rubber (mushroom type) ear tips
- Adult and pediatric diaphragms



CAT#	DESCRIPTION	BOX
6-0028-00	Sprague Rappaport, Black	1
6-0028-50	Sprague Rappaport, Black	10



# Accessories



## Monascope / Duoscope

### ComfortFIT Monascope

- Disposable
- Latex free

The transmission of sound through this system (earpiece, tubing, and stethoscope head) is dependent on a complete seal in the system. At the patient's end this is not difficult to accomplish. At the anesthetist's end this is more difficult due to the irregular and variable shape of the ear canal.

The Comfort FIT earpiece answers this need perfectly. The foam is compressed, inserted into the ear canal and slowly expands to make a perfect, low pressure seal with the canal wall. This seal is critical due to the low volume of the sounds transmitted. The low pressure against the skin in the ear canal makes the earpiece more comfortable which encourages more vigilant use. Tubing is four feet in length.

Connects to:

- Esophageal Stethoscope
- Wenger Chest Piece




CAT#	DESCRIPTION	BOX
6-1030-25	Comfort FIT™ Monascope	25

## DUOSCOPE

- Latex free

A standard binaural stethoscope piece designed specifically for the anesthesiologist to provide maximum comfort and efficiency. Tubing from Y-piece to female luer connector is extra-long (36"), permitting considerable latitude of movement around patient without detachment.

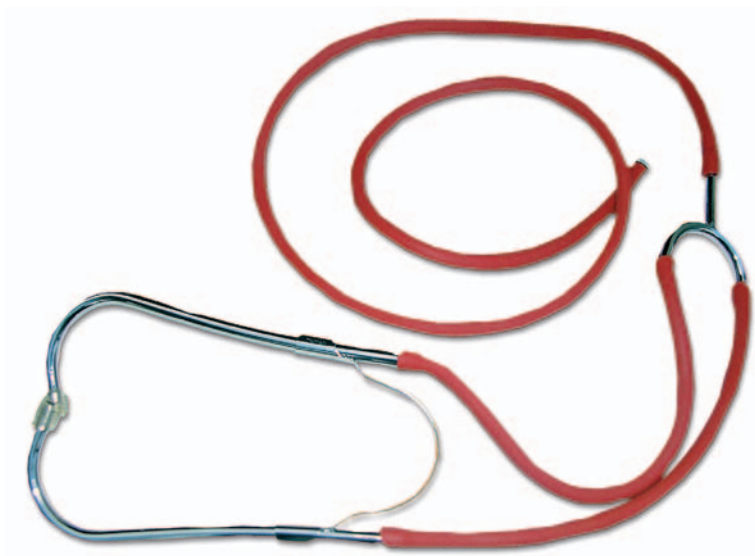

CAT#	DESCRIPTION	BOX
6-0030-00	Duoscope	1



# Accessories



## Chest Pieces (Precordial Stethoscopes)

### CHEST PIECES

#### INFANT CHEST PIECES

Lightweight, stainless steel infant chest pieces feature a hollow design and rounded, smooth chest surface, making them ideal for use with newborns and infants.

#### WENGER CHEST PIECES

Solid chest pieces made from stainless steel. They have relatively deep sound chambers and are supplied with male luer connectors for attachment to ear pieces.


- All chest pieces have a luer slip male connector
- Stainless steel - no chrome plating to flake or peel



CAT#	DESCRIPTION	BOX
6-0025-00	Newborn Chest Piece, 7/8" diameter	1
6-0025-01	Infant Chest Piece, 1 1/8" diameter	1
6-0025-02	Child Wenger Chest Piece, 1 1/8" diameter	1
6-0025-03	Adult Wenger Chest Piece, 1 7/16" diameter	1

### DOUBLE ADHESIVE DISK

CAT#	DESCRIPTION	ROLL
6-0025-25	Double Adhesive Disk	1

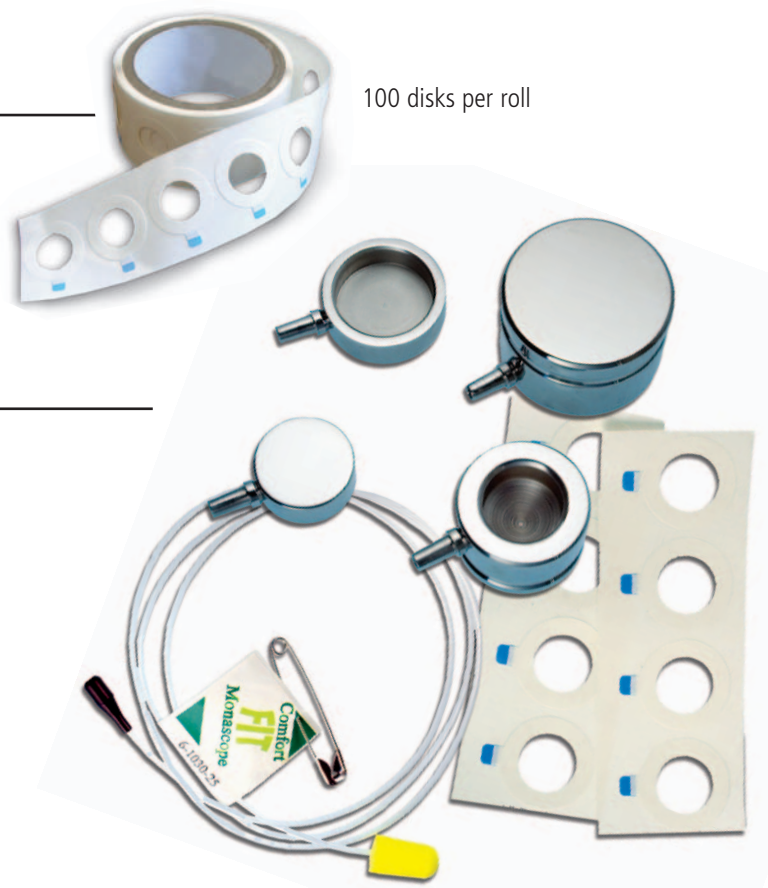

100 disks per roll

### CHEST PIECE KIT

KIT INCLUDES:

- 1 each adult, child, infant, & newborn chest piece
- 8 Double adhesive disks
- 1 Comfort FIT™ disposable monoscope

CAT#	DESCRIPTION	KIT
6-0025-05	Chest Piece Kit	1



# Accessories



## Eye and Tooth Protector / Nose Clips

### IGuard™ EYE PROTECTOR

During general anesthesia, protect your patient's eyes against corneal abrasion, lacerations and scleral hemorrhage with the IGuard eye covers. These rigid, clear, plastic protective covers are supplied with a non-allergic self-adhesive foam cushion for fast and accurate application.

- Single patient use
- Sterile
- Latex free

CAT#	DESCRIPTION	BOX
9-0210-00	IGuard Eye Protector	25



### DISPOSABLE TOOTH PROTECTOR

Helps to prevent broken or chipped teeth during intubation, bronchoscopy, laryngoscopy, tonsillectomy, and tooth extraction. Minimizes tooth injury: One third of all anesthesia malpractice suits involve teeth damaged from intubation.

- Adult size
- Disposable
- Latex free

CAT#	DESCRIPTION	BOX
9-3041-03	Adult Tooth Protectors	50



### DISPOSABLE NOSE CLIP

- Latex free
- Disposable

CAT#	DESCRIPTION	BOX
9-0134-44	Nose Clips	100

